


مؤسسة تخصص العناية للتجارة CARE SPECIALITY EST.



0 Ζ Ζ



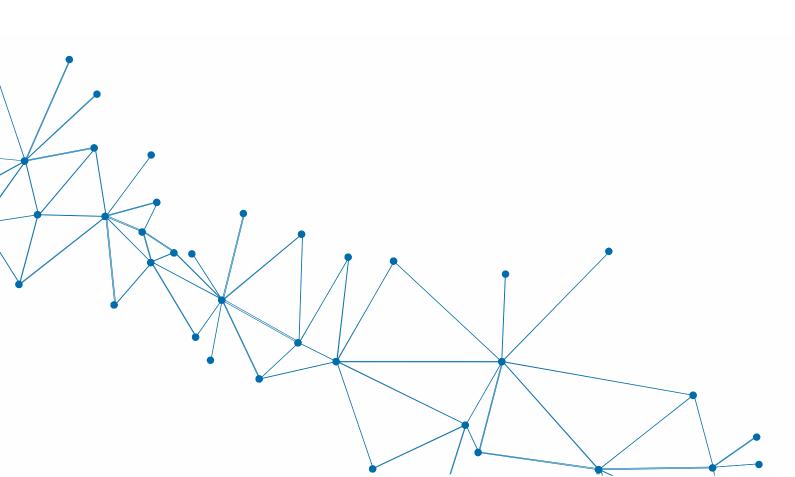





your direction

DENTAL LIGHT

WASHBASH

STERUXER

GIVE US CHANCE AND WE WILL PROVE OUR EFFICIENCY

1-RAY

DENTALMASY

DENTAL GLASSES

DENTAL INSTRUMENTS

K-RAY MACHINE

TOOTHIMPLANT

TOOTHBRUSH

DENTIST CHAIR



Cordless Endo MOTOR

- The contra angle can be rotated in 340 degrees, making it easy to watch the LCD during the treatment. Quick and easy programmable modes make easier to operate.
- Miniature Head: Provides excellent visibility while working.

- Quick and simple programmable modes for torque, speed, continuous rotation/ reciprocation modes

- Auto Torque Reverse: Reduces the possibility of instrument separation inside the canal. - Reciprocation mode for regular rotary files like Neoendo Flex.
- Adjustable for both left and right hand converting, self-calibration, volume adjustment etc.
- Rotational Speed Range -120-650 RPM
- Torque Range –0.5 4.0 N cm
- Battery Capacity –1600 m Ah
- Auto Start-stop- Motor starts moving when the file enters the canal and stops when withdrawn.
- Apical Reverse- When the file reaches pre-set apical reverse point.
- Apical Slow Down- To preserve the cement-dentinal junction and avoid file-separation, speed of motor slows down when the file is near apical region
- Apical Torque Reduction- In addition to apical slow down, torque reduction at apical third prevents chances of file separation.



APEX LOCATOR

- Compact size and stable design allow for excellent flexibility during procedures.
- Based on advanced multiple frequency network impedance measurement technology and automatic calibratina ensures the measurements are accurate.
- OLED display, brighter and cleaner surface.
- Battery Capacity ; 1600 m Ah
- Screen ; 3.5" LCD



Endo motor with built-in **Apex locator**

- Reciprocation Adjustable speed and Adjustable angles from 30 to 370 degrees.
- Built-in Apex Locator: integrate with motor function or use as standalone apex locator.
- Ambidextrous Display Screen can be rotated for left or right-hand use. Using OLED screen.
- Pre-loaded with all major file systems, which can be further customized.
- Auto Start and Stop Motor starts moving when the file enters inside the root canal, and stops when withdrawn.
- Apical Reverse: When the file reaches pre-set apical reverse point, user can define whether the file starts moving in reverse direction.

- Apical Slow Down: to preserve the cement-dentinal junction and avoid file separation, speed of motor slows down when file is in apical region.

Technical

- Rotational speed range: 120 to 1000 rpm
- Torque range: 0.5-4.0 N/cm
- Battery capacity: 1500 m/Ah



LED Curing Light

- 3 Seconds for Curing
- Light Intensity: 1600 mW/cm2 High Intensity, with Three working modes: LOW, MIDDLE, HIGH.
- Stable output of light intensity, the solidification effect is not affected by the consumption of remaining power.
- Wireless. Full charge used for more than 400 times continuously under 10 sec working time mode.

- Easy-to-change battery. Lithium battery with CE certification, more reliable, no explosion. Battery on the market without sound charge management will easily cause explosion under high temperature with unguaranteed quality.

- Low standby power consumption with 40 days' standby time.
- Ergonomics design, more convenient.
- Unique light hood holder design, high efficiency light filter.
- Working model: Full, Ramping, Pulse.
- Time setting: 5s, 10s, 15s, 20s
- Light hood

Studies prove that Blue ray will cause macular degeration to retina and eyes.

Light hood is mode of selected quality materials, preventing eyes from blue ray injury.

Digital Metal Ultrasonic Cleaner

Process works by sending ultrasonic sound waves through water, thus creating millions of tiny bubbles that literally lift the dirt and grime off your dental device. Features:

- -External drainage system makes easy change of rinse liquid.
- -Low noise design makes quite operation.
- -Total high precision plastic cover
- -Rinse time can be set from 1 to 60 min
- -The even temperature from 20-80
- -Power is automatically cut off when is lack of water.
- -Low noise design makes quiet operation.
- -Specialized basket is designed for delicate minute devices that should be free from
- friction and braise, for example, high-speed and low-speed hand piece turbines.
- -Clear operation plate with LED display makes handy operation.

Technological specifications: -Volume: 5L -Temperature range: 20-80°C -Time range: 1-60 minutes

SEALING MACHINE

- PTC Heating bar, automatic thermostat, long life
- Automatically alarm when wrong operation or spare parts broken
- High quality plastic shell design, generous and practical
- Easy to operate, long time, the power is automatically turned off

Feature:

- Stable function
- Low rate of repair
- The build-in double blades are safe to operate with smooth incision
- Industrial-grade double-open blade, safe and sharp, long life for using



Whitening Accelerator

The best laser treatment products available on the marketplace. In case you have apparently observable and tough stains from food, drinks or smokes, then is particularly effective. Even if you have fluoride or age-related stains, will bring back your teeth to as naturally whitened as they may be. Teeth-whitening not just makes the teeth brighter, but also removes discolorations which are due to several factors. Diet, medications, coffee, and cigarettes - they all affect the enamel of the teeth.

The Whitening Accelerator system is completely safe and allows teeth whitening also in patients with sensitive teeth. The last effect depends mainly on the nutrient composition of the tooth and stabilizes about within a week following the procedure. Following this time, it is likely to replace fillings, crowns, and other prosthetic structures which do not lighten with whitening and are found in the teeth that show when you smile.

This filtration system eliminates harmful heat and UV light in the halogen light source and brings the complete power of the light into the front of the lamp. The 12 blue LED lights unify using the halogen light source in the front of the lamp head providing the most innovative blue light in a wavelength of 480-520 nanometers.

In the middle of the Advanced Light Technology is ultrasonic technology that's used to assist in the oxidation procedure. Operates at a lower temperature compared to similar products using the latest optical technology thus preventing irritation of the nerves of the teeth.

When used in combination with our proprietary whitening formula which includes Hydrogen Peroxide (H2O2) and Silicon Dioxide (SiO2) that is a catalyst, along with other special ingredients, the end result is not just whiter teeth but shinier teeth with reduced prep and seat time.

Features:

- Halogen-powered teeth whitening accelerator, LED curing light, and Low-Level Laser Therapy (LLLT) option
- Powerful halogen light output, clinically proven to provide superior results to LED and plasma arc systems
- Unique Lightbridge light technology
- Dual-arch whitening features gentle, normal, and intense light-output settings for the perfect treatment setting with each individual patient
- Single-tooth whitening option included for precision tailored treatment results
- Three (3) Minute treatment complete warning timer.
- Custom fitted cover for easy storage.
- 97% clinical rating on dentaladvisor.com
- wireless remote and patient call button
- advanced air ionizer to remove heat and Ultraviolet output
- curing light has three settings to match the needs of your customers
- motorized arm

Unprecedented Innovation and Flexibility in Whitening Technology:

- Upgraded light filtration system containing over 200,000 optical fibers, to filter out all harmful heat and ultraviolet (UV) light.
- Dual-arch power whitening option can be set to High, Medium, or Low light output settings for a tailored treatment that best suits each individual patient.
- Single-tooth whitening option included for precise, tailored treatment results.
- Smart Remote Control that allows patient awareness of procedure time.
- Built-in purifier for a cleaner and healthier environment.

Innovative and Power-Packed System:

- The improved motorized arm that allows for yet easier adjustment of the accelerator head for a more comfortable whitening procedure.
- State of the art, full colored touch screen with a newly added multilingual feature.
- Built-in instruction mode, that allows for a customized whitening procedure based on the client's needs.
- Built-in, powerful LED curing light, featuring a dual wavelength output to cure any dental composite on the market.
- Radial arm with 360 range of motion.
- Motorized height adjustment can be controlled from the TUI control panel and accelerator head.
- Advanced temperature cooling and monitoring system designed to prevent a system from overheating.



Food-Grade Material

All the parts that contact with liquid are made of 304 stainless steel, including inner tank, inner cap, condensing coil, water outlet, and water inlet. Featuring BPA-free plastic container and Aluminum cooling fan.





Convenient Button

A simple button easily realizes start function without having to plug or unplug, featuring an indicator light to show the status. Dry burning-resistant protection function included.

Water Distiller





Considerate & Portable Design

The upgraded filter is made of 304 stainless steel; a unique and portable handle on the top for easy carry; a built-in thermostat to prevent overheating and ensure auto-off (when the temp is over 320°F).



High Efficiency

With high efficiency, it takes only about 2-3.3 hours to make 4L (1.06 US gal) distilled water. 50Hz 750W heating element makes distillation speed over a quarter gallon per hour, distills up to 12 gallons per day.

Effective Heat Dissipation

The plastic wind leaf change to aluminum wind leaf, improving the cooling effect. The discharge cap changes from plastic to 304 stainless steel, improving the heat dissipation further.





COOLING FAN

AIR VENT

Dental Chairs S

Deluxe removable cleaning ceramic spittoon Button protective plastic cover with a three bit memory a electric dental chair Tunable high-grade LED operating light light Auto spitton flush and cup filler control system high-grade articulated headrest 24V noiseless Denmark DC motor Air brake instrument try Saliva ejector with high and low suction 3-way syringe Anti rust body LED x-film viewer Dentist Stool A Full function foot control





Dental Chairs P

Main Features:

- Luxury shadow-less LED light 20,000Lx with Color temperature 6500K
- Spitton cabinet can move laterally 45 degree with Luxury 2-Colors
- 9-memories, working status LCD display
- Full control from main tray (Positions, Memories, Cup/Spitton water, LED light, LED X-Ray film viewer)
- Rotatable handpiece holder (with 4-modules for handpieces + Air/Water Syringe module)
- 24V DC Denmark motor
- Multi-function foot control with Luxury metal design
- Strong, weak a saliva suction system
- Doctor Stool: Large type Natural Leather
- LED X-Ray film viewer
- Built-in tissue box
- Natural Leather with fantastic colors
- Cup holder
- Glass tray

Intraoral X-Ray Sensor

- World's thinnest and most economical intraoral X-ray sensor ever- 4.4mm.
- AED Technology (Automatic Exposure Detection) -from chest radiography
- Compatible with almost all DC and AC X-Ray generators
- Digital images come automatically when your X-Ray is switched on
- APS CMOS Technology -from aerospace industry
- High signal-to-noise ratio with maximum details
- Long life and low power consumption
- Direct Deposition CSI Technology -from chest radiography
- Low dose radiation for safety
- Less exposure time but higher resolution
- All in one the most economical X-Ray sensor in the world

Unique designs

- Cupped cable connection design brings easy-position for finger
- Robust Cable Connection70000+ Bending RA test passed
- User friendly 4 step software
- True resolution of 20 lp/mm
- AD conversion- 16 bit
- IP68- highest water and dust protection rating
- Tablet with pre-loaded software
- Interest free EMI
- 1 / 10 Cost per patient compares to film
- 3s to complete whole imaging process
- 50000+ Exposure with no quality loss

Ultra-thin body design Available with and without tablet



Unique designs

Ergonomics design Cupped cable connection design brings easy-position for finger

Robust Cable Connection

Ultra-thin body design



Portable Dental X-Ray

Product Features:

The device is DC high frequency portable dental x-ray unit, small, light and nearly no radiation.

There are manual buttons installed on the surface of shell, operate easily.

The unit mainly suitable for oral pre-treatment diagnosis of internal

organization structure and root depth and so on, is indispensable in daily life

clinic, especially for dental implant surgery.

The battery is durability, it can take about 500 images after full charged, and it can be charged and discharged for 1000 times in its life.

It can be connected with intra-oral digital x-ray imaging system.

Specifications

Tube voltage:60KVTube current:1.5mAExposure time:0.1-2sFrequency:30KHzTube focus:0.3mm×0.3mmTube focus:0.3mm×0.3mmBattery: DC14.8V 6400mAhCharger input voltage:AC100V-240V±10%Ambient temperature:5-40°CRelative humidity:<75%RH</td>

Net Weight:2.5kg

dental disposable





dappen dish	CS - 1122 CS - 1	CS - 1123
impression plastic tray	CS - 1125	CS - 1126
CS - 1127 high volume suction	CS - 1128	CS - 1129
CS - 1130	CS - 1131	CS - 1132 alcohol swabs Alcohol swabs Provincement Pro











6	
	$\left[\mathcal{O} \right]$



